


**APPROVED IN
EUROPE AND CANADA**

BackStop™ for kidney stone removal


A new, completely atraumatic device to prevent stone migration and stone fragment retropulsion during ureteroscopic lithotripsy

- Prevents retrograde stone migration
- Prevents scattering of all stone fragments
- Does not interfere with the lithotripter
- Easily used in ureteroscopy procedure
- Available in clear or blue


BackStop™ is delivered above the stone.

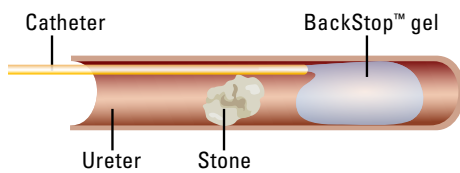


The catheter is withdrawn to free the operating field.

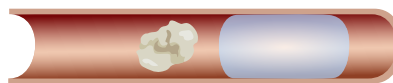

BackStop™ prevents stone migration during lithotripsy.

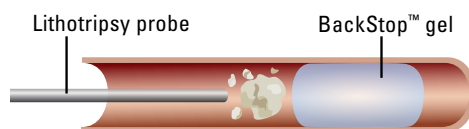

BackStop™ is removed via conventional irrigation.

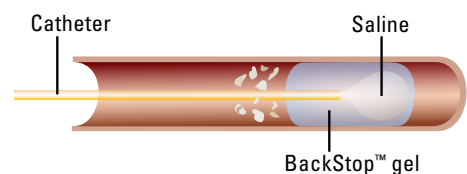

Designed and approved to prevent retrograde stone migration during ureteroscopic lithotripsy.



BackStop™ with the BackStop™ Injector and Catheter

Migration of stones and stone fragments during ureteroscopic lithotripsy can lengthen and complicate the procedure, or necessitate subsequent procedures. Traditional baskets or wires may cause trauma to the ureter, interfere with the lithotripter and allow migration of smaller stone fragments.

BackStop™ is a water soluble gel that is delivered above the stone in the ureter and acts to prevent retrograde stone migration during ureteroscopic lithotripsy. Once delivered, the catheter is removed, thus freeing the operating field. BackStop™ is completely atraumatic and will conform to any shape ureter. It has no mechanical parts and will not interfere with the lithotripter. BackStop™ is compatible with all forms of ureteroscopic lithotripsy, including laser and mechanical lithotripters.

When lithotripsy is complete, BackStop™ is dissolved via conventional irrigation and the ureteral occlusion is reversed.

Ordering Information - contact europe@pluromed.com or canada@pluromed.com

catalog numbers		
BackStop Syringe	BackStop Injector	BackStop Catheter
Catalog #60-0250: 2.5 mL BackStop™ Syringe	Catalog #60-K25: BackStop™ Injector – 2.5mL	Catalog #60-0000-3100: 3F BackStop™ Catheter
Catalog #60-0500: 5.0 mL BackStop™ Syringe -or- Catalog #65-0500: 5.0 mL BackStop™ Blue Syringe	Catalog #60-K50: BackStop™ Injector – 5.0mL	Catalog #60-0000-3100: 3F BackStop™ Catheter -or- Catalog #60-0000-5100: 5F BackStop™ Catheter